- EN Quick Guide
- FI Pikaopas
- SV Snabbguide
- DA Quick Guide
- NO Hurtigguide
- EL Γρήγορος οδηγός

iCare PATIENT2
iCare EXPORT

Quick Guide

icare

iCare PATIENT2 & EXPORT
Quick Guide for
mobile phones and PC

Read this Quick Guide before using iCare PATIENT2. Follow the instructions step by step.

If you need help in using the PATIENT2 app, contact your healthcare professional or iCare distributor.


English

PATIENT2 Quick Guide for Android

Step		Page
1	Start using the iCare PATIENT2 app	5
2	Transfer measurements and Log in	6
3	Graph functionality	8

PATIENT2 Quick Guide for iOS

Step		Page
1	Start using the iCare PATIENT2 app	11
2	Pair the tonometer with your mobile device	12
3	Transfer measurements and Log in	14
4	Graph functionality	16
5	How to switch on Bluetooth on the Tonometer	18

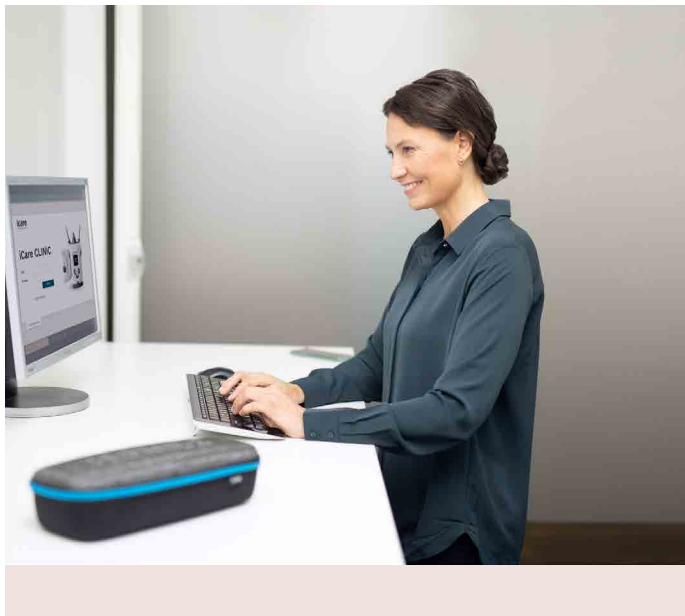PATIENT2 Quick Guide for Android

Patients can use the iCare PATIENT2 for viewing their IOP measurement results and for transferring results to the iCare CLINIC or to the iCare CLOUD. Healthcare professionals can use the iCare PATIENT2 app only for transferring IOP measurement results to the iCare CLINIC.


Make sure that the mobile device with the iCare PATIENT2 app installed is set to the correct time and date.


Before taking measurements, update the tonometer's time to your local time manually from the tonometer's settings or automatically by connecting the tonometer to the iCare PATIENT2 application or to the iCare EXPORT software.


2

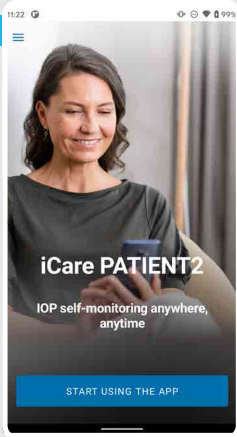3

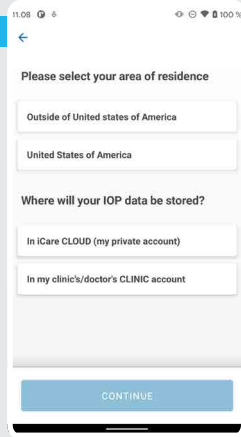4

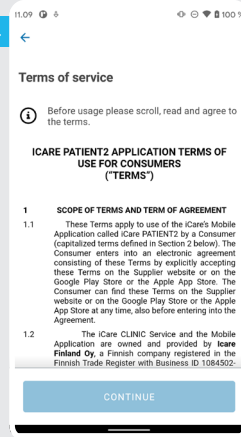5

1

Download the iCare PATIENT2 app from Google Play Store.

2

Press **START USING THE APP** button.

3




Select your area of residence. Select where your IOP data will be stored.

4


Scroll down and read through the Terms of service and agree to them.

5

You are now ready to start using your iCare PATIENT2 app.

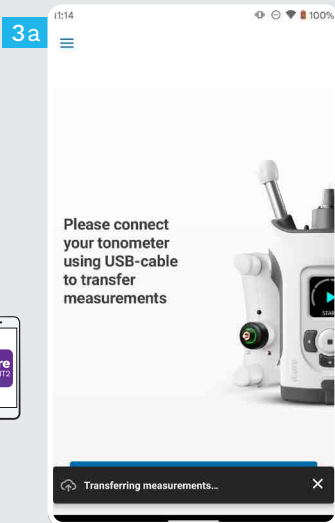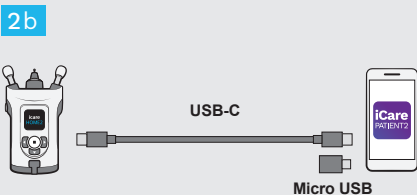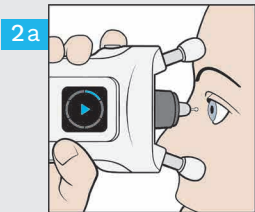
- 1 Make sure your mobile device's screen and keys are not locked.
- 2 After taking measurements, connect your tonometer to your phone using the USB cable provided. You can use the provided adaptor if needed.
- 3 Once the USB cable is connected, the measurement data is transferred automatically from the tonometer to the cloud. Do not disconnect the cable while the measurement data is being transferred.

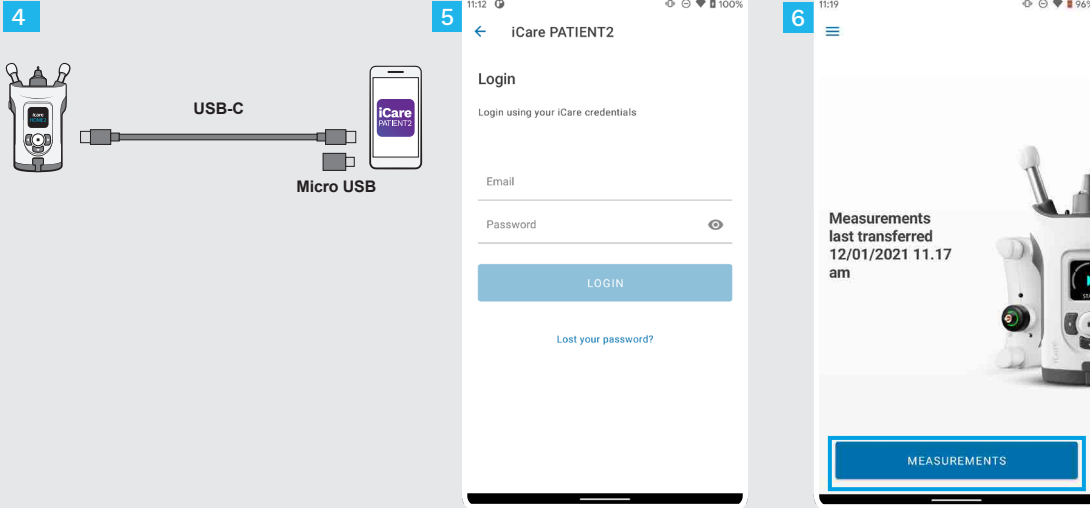4 Disconnect the USB cable after the measurement transfer is completed.

5 Log in with your iCare credentials. For information on how to obtain iCare credentials, read the Getting Started Guide in the tonometer's carrying case.

6 Tap the **MEASUREMENTS** button to see your IOP measurement results.

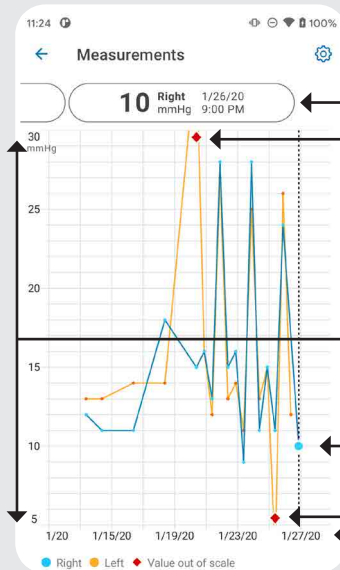Settings 


Selected measurement

Red indicator  shows the measurement values outside of the graph scale.

Graph scale can be adjusted from .

Selected measurement

Measurement date and time

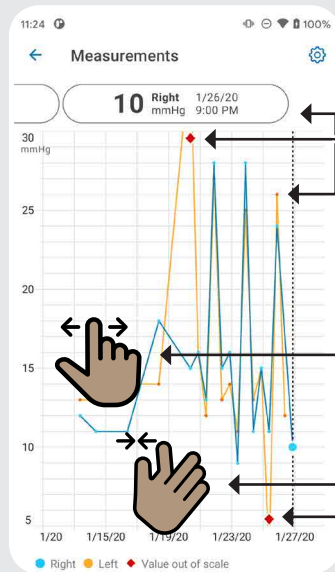
To see the result of each measurement, scroll the pill menu at the top of the screen or tap the dot on the graph.

Tap **◆** to see the measurement value that is outside of the graph scale.

Scroll on the graph or on the pill menu to see your previous or next measurements.

Pinch on the touch screen with two fingers to zoom in to see fewer IOP results on the graph. Zoom out to see the whole graph.

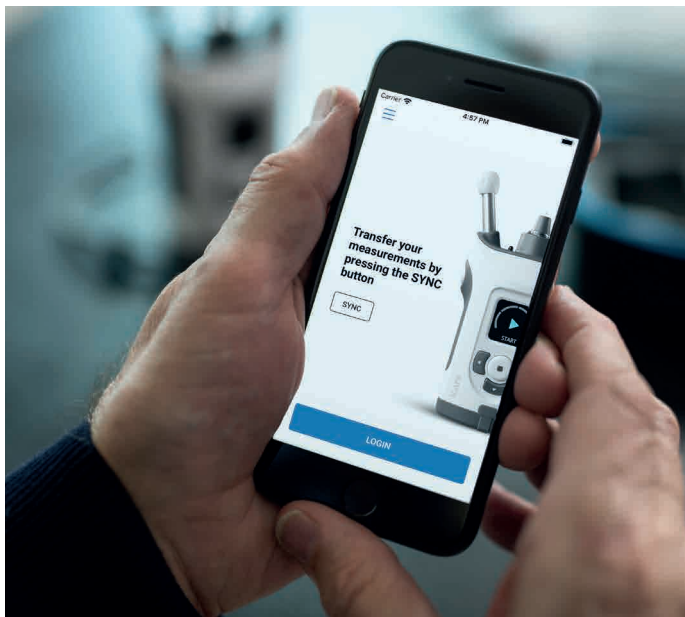PATIENT2 Quick Guide for iOS




Patients can use the iCare PATIENT2 for viewing their IOP measurement results and for transferring results to the iCare CLINIC or to the iCare CLOUD. Healthcare professionals can use the iCare PATIENT2 app only for transferring IOP measurement results to the iCare CLINIC.

Make sure that the mobile device with the iCare PATIENT2 app installed is set to the correct time and date.

Before taking measurements, update the tonometer's time to your local time manually from the tonometer's settings or automatically by connecting the tonometer to the iCare PATIENT2 application or to the iCare EXPORT software.

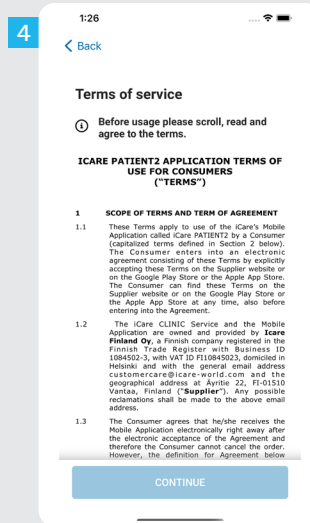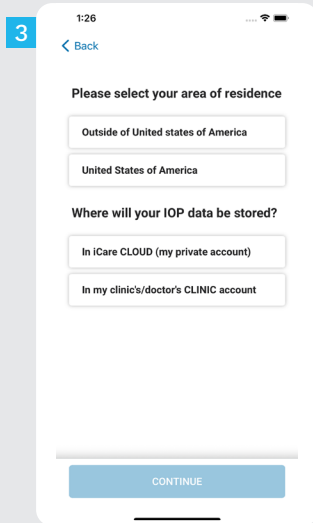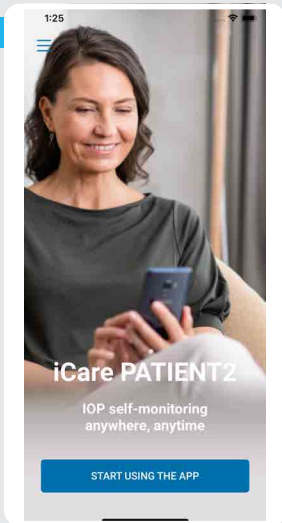1 Download the iCare PATIENT2 app from App Store.


2 Press "START USING THE APP" button.

3 Select your area of residence. Select where your IOP data will be stored.


4

Scroll down and read through the Terms of service and agree to them.


3a

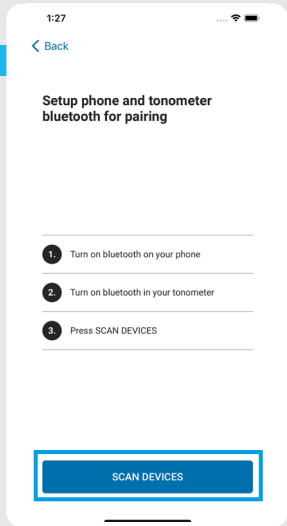3b

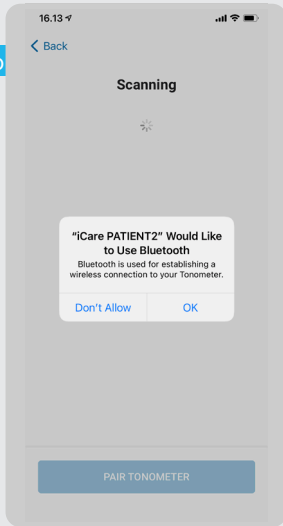4

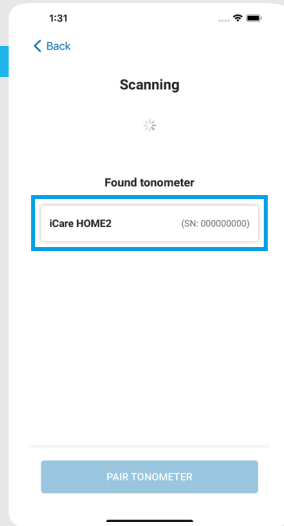1

Check if Bluetooth is on on your tonometer and if not, turn it on. See Chapter 5 for instructions.

2

Go to your smart phone's settings and switch Bluetooth on.

3




Tap **SCAN DEVICES**. If the app asks for a permission to use Bluetooth, tap OK.

4

Tap on the tonometer with the serial number.

5



Tap **PAIR TONOMETER**.

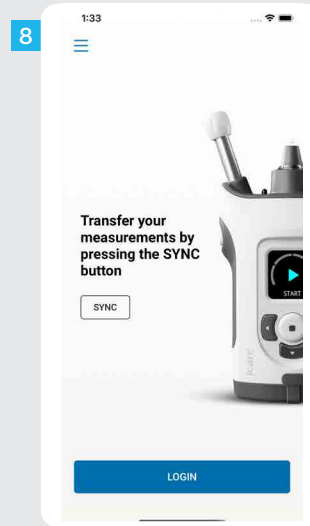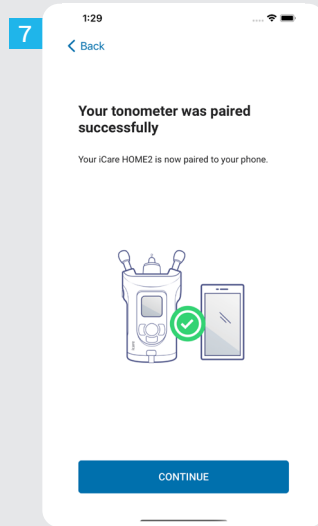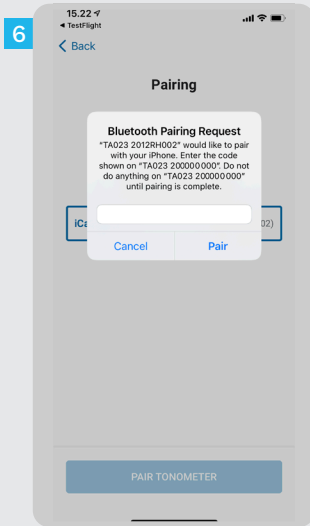
6 Enter the PIN code shown on the tonometer's screen within 30 seconds and tap **Pair**. Otherwise, you need to restart the pairing process.

7 When the tonometer has been successfully paired with your mobile device, tap **CONTINUE**.

8 You are now ready to start using your iCare PATIENT2 app.

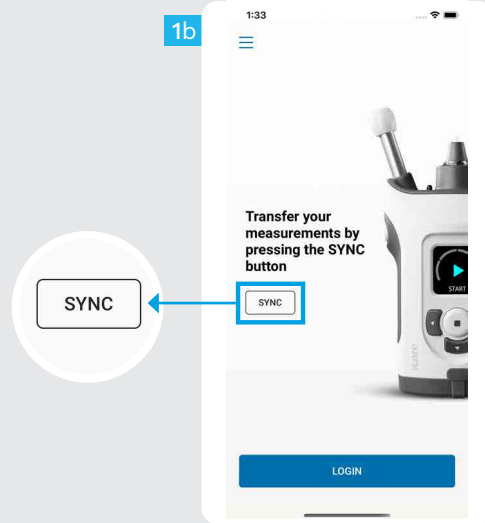
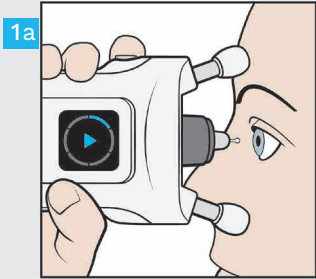
1 After measuring your IOP you can use the **SYNC** button to transfer the measurement results to the cloud.

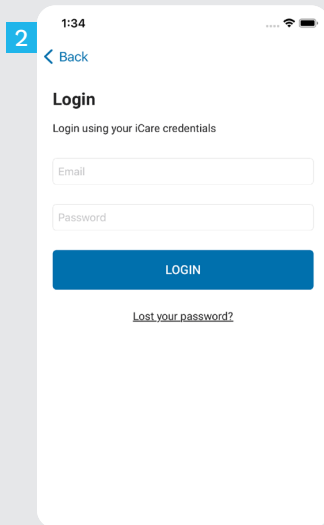2 Log in with your iCare credentials. For information on how to obtain iCare credentials, read the Getting Started Guide in the tonometer's carrying case.

3 After logging in, you can see your IOP measurement history by tapping the **MEASUREMENTS** button.

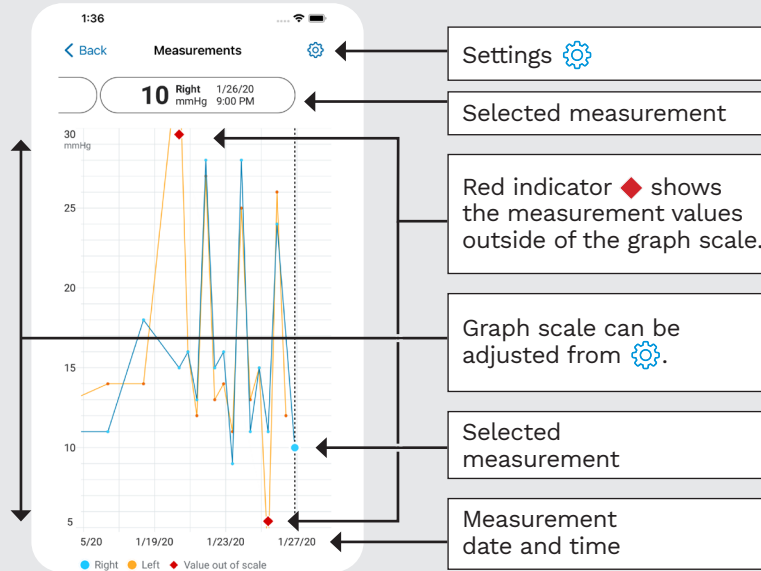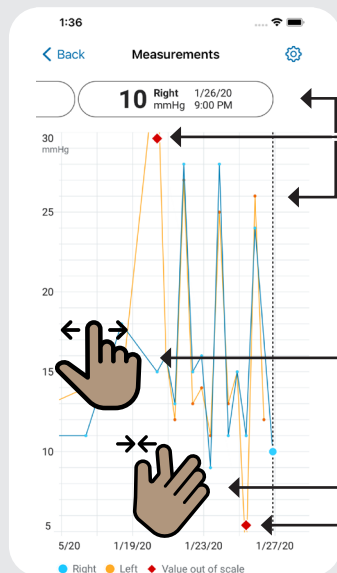To see the result of each measurement, scroll the pill menu at the top of the screen or tap the dot on the graph.

Tap  to see the measurement value that is outside of the graph scale.

Scroll on the graph or on the pill menu to see your previous or next measurements.

Pinch on the touch screen with two fingers to zoom in to see fewer IOP results on the graph. Zoom out to see the whole graph.

Hold the tonometer vertically for these instructions.

The Bluetooth on your tonometer is on by default as a factory setting. In case you want to check if the Bluetooth is on on the tonometer or turn it on, follow the instructions below.

- 1 Press  to switch the tonometer on.
- 2 Press  until you see **SETTINGS** on the tonometer screen, and press .
- 3 If Bluetooth is off, Press  to enter the **BLUETOOTH** menu.
- 4 Press  until you see **TURN ON**, and press .

iCare EXPORT


Quick Guide

icare

icare EXPORT

Quick Guide (For Windows only)

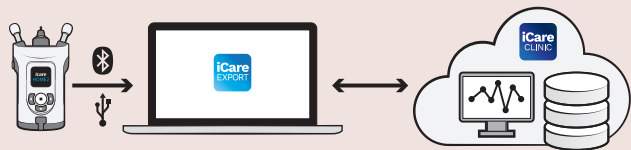Read this Quick Guide before using iCare EXPORT. Follow the instructions step by step.

If you need help in using the EXPORT software, contact your healthcare professional or iCare distributor.

EXPORT Quick Guide

Step

Page

- | | | |
|----------|---|----|
| 1 | Download and Set up iCare EXPORT | 21 |
| 2 | Connect and send measurements using USB cable | 23 |
| 3 | Connect and send measurements via Bluetooth | 24 |



Patients can store their IOP measurement data to a private account in the iCare CLOUD, if their tonometer is not registered to a healthcare professional's iCare CLINIC account. For information on how to set up a private account, read the Getting Started Guide in the tonometer's carrying case.

Make sure that the computer time and date is correct when using with iCare EXPORT.

Before taking measurements, update the tonometer's time to your local time manually from the tonometer's settings or automatically by connecting the tonometer to the iCare PATIENT2 application or to the iCare EXPORT software.

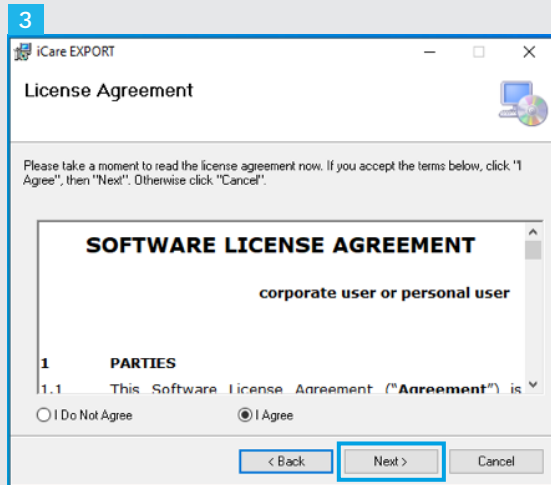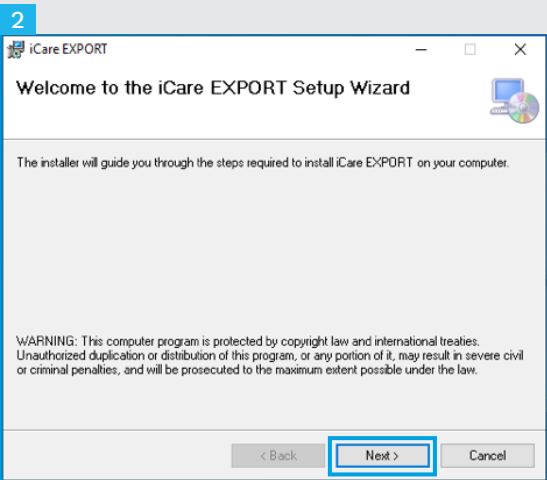1 Download iCare EXPORT from the following website:

<https://www.icare-world.com/apps/export/>.

Double click on the iCare EXPORT Setup.exe file.

2 Read the instructions on the Setup Wizard and click "Next".

3 Accept the software license agreement by selecting "I Agree". To continue, click "Next" until the installation is completed.


The image displays three sequential screenshots of the iCare EXPORT software installation and interface:

- Step 4:** A screenshot of the "iCare EXPORT" window titled "Installation Complete". The text inside reads: "iCare EXPORT has been successfully installed. Click 'Close' to exit." At the bottom, there are three buttons: "< Back", "Close", and "Cancel".
- Step 5:** A screenshot of the desktop showing the "iCare EXPORT" icon, which is a blue square with a white lowercase 'i' and the text "iCare EXPORT" below it.
- Step 6:** A screenshot of the iCare EXPORT application window. The title bar says "iCare EXPORT". The main content area says "Connect device to continue" and "Please connect your tonometer to computer with USB-cable or Bluetooth". Below this is a table with columns: "Measurement time", "IOP (Right eye)", "Quality (Right eye)", "IOP (Left eye)", and "Quality (Left eye)". The table is currently empty. At the bottom left, there is a "Select BT Device" dropdown menu showing "No device" and a "Copy" icon. At the bottom right, there are buttons for "Save local copy", "Open local copy", and "Clear table". A language selection menu is open on the right side, listing various languages including Danish (DK), Deutsch (DE), English (GB), English (US), Español (ES), Suomi (FI), Français (FR), Italiano (IT), 日本語 (JP), 대한민국 (KR), Norsk (NO), Polski (PL), Português (BR), Pycckî (RU), Slovenčina (SI), Svenska (SE), 繁體 (TW), and English (GB).

4 When installation is complete, click “Close”.

5 Double click the “iCare EXPORT” icon on the desktop or find the program under “Windows Start” menu, under “Icare” folder and start the program.

6 You can select the language from the right bottom menu.

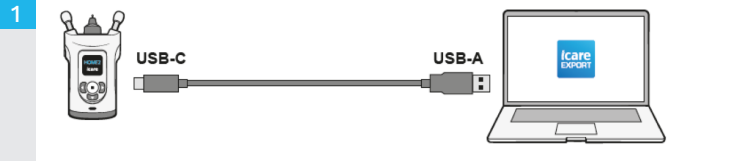2

icare EXPORT

LOG IN →

Measurements saved!

All measurements were successfully posted (Organization 1 in London Demo). The measurements were cleared from the device.

Measurement time	IOP (Right eye)	Quality (Right eye)	IOP (Left eye)	Quality (Left eye)
1/7/2021 10:04:17 AM			16	GOOD
1/7/2021 10:04:06 AM	14	EXCELLENT		
1/7/2021 10:03:59 AM			14	GOOD
1/7/2021 10:03:52 AM	14	GOOD		

Select BT Device

2028RH005 S/N



Copy

Save local copy Open local copy Clear table

1 Connect your tonometer to your PC using the USB cable provided.

2 The tonometer sends the measurement results to the server automatically.

NOTE! that once you have transferred the data, the measurement results are automatically removed from the tonometer's memory.

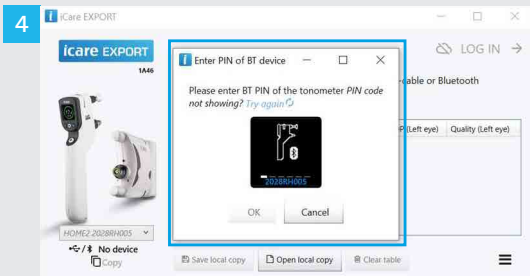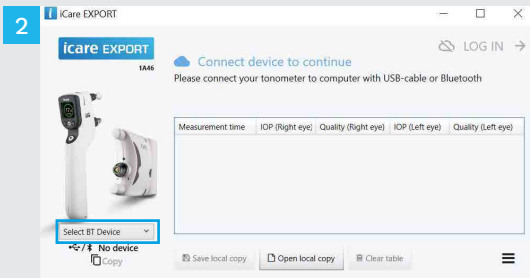In iCare EXPORT application

For these instructions, make sure that the tonometer's and the PC's Bluetooth is on. For instructions, see chapter 5 How to switch on Bluetooth on the Tonometer.

- 1 Make sure that the tonometer is on before you connect the tonometer to iCare EXPORT.
- 2 Click "Select BT Device" on the left bottom.
- 3 Select your tonometer's serial number from the drop-down menu. You can find the serial number (SN) on the label on the back of your tonometer.
- 4 Enter the PIN code that will appear on your tonometer's screen and press "OK"

1 **icare EXPORT**

icare EXPORT 1A46 LOG IN →

Measurements saved!
All measurements were successfully posted ([Organization 1 in London Demo](#)) The measurements were cleared from the device.

Measurement time	IOP (Right eye)	Quality (Right eye)	IOP (Left eye)	Quality (Left eye)
1/7/2021 10:04:17 AM			16	GOOD
1/7/2021 10:04:06 AM	14	EXCELLENT		
1/7/2021 10:03:59 AM			14	GOOD
1/7/2021 10:03:52 AM	14	GOOD		

Select BT Device
2028RH005 S/N Copy

Save local copy Open local copy Clear table

1 The tonometer connects to the iCare EXPORT application and sends the measurement results to the iCare CLOUD or your healthcare professional's iCare CLINIC account automatically.

NOTE! *that once you have transferred the data, the measurement results are automatically removed from the tonometer's memory.*

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Äyritie 22
01510 Vantaa, Finland
Ph. +358 9 8775 1150
info@icare-world.com

Icare USA, Inc.

4700 Falls of Neuse Rd. Ste 245
Raleigh, NC. 27609
Ph. +1 888.422.7313 | Fax +1 877.477.5485
infoUSA@icare-world.com

For further information,
read the instruction manual.

www.icare-world.com

icare

iCare is a registered trademark of Icare Finland Oy, Centervue S.p.A., Icare Finland Oy, and Icare USA Inc. are parts of Revenio Group and represent the brand iCare. Not all products, services or offers referenced in this brochure are approved or offered in every market and approved labeling and instructions may vary from one country to another. Product specifications are subject to change in design and scope of delivery and as a result of ongoing technical developments.